

PRO3-P™ and PRO3-C™

Amniotic Membranes



DONATED HUMAN TISSUE

RESTRICTED TO USE BY OR ON THE ORDER OF A LICENSED HEALTHCARE PROVIDER (physician, dentist, podiatrist, optometrist, nurse practitioner or physician assistant). NOT INTENDED FOR VETERINARY USE.

80-129 Rev. 07
P01-IFU-0004 Rev H

DESCRIPTION

PRO3 amniotic membrane (a product family consisting of PRO3-P™ amniotic membrane and PRO3-C™ amniotic membrane) is a semi-transparent, collagenous membrane obtained with consent from mothers during childbirth. PRO3 amniotic membrane is derived from placental tissue. PRO3 amniotic membrane is processed using aseptic techniques and dehydrated. The allograft is aseptically packaged in a tear pouch within a peel pouch and provided sterile.

INTENDED USE

PRO3 amniotic membrane is intended for use as a tissue barrier or wound covering.

CONTRAINDICATIONS

PRO3 amniotic membrane should not be used in areas with severe vascular compromise or with active or latent infection or in a patient with a disorder that would create an unacceptable risk of post-operative complications.

DONOR ELIGIBILITY

The tissue that comprises PRO3 amniotic membrane is acquired from a qualified donor and processed using aseptic techniques in accordance with federal, state, and/or international regulations and to the standards of the American Association of Tissue Banks. The donor is screened and tested for communicable disease risks and other exclusionary medical conditions. The results of the donor screening and testing have been reviewed by the Medical Director of Vivex Biologics, Inc., and the donor has been deemed suitable for transplantation.

Communicable disease testing is performed by an FDA-registered laboratory certified to perform such testing on human specimens under the Clinical Laboratory Improvement Amendments of 1988 (42 U.S.C. 263a) and 42 CFR Part 493, or that has met equivalent requirements as determined by the Centers for Medicare and Medicaid Services in accordance with those provisions. Results from the following infectious disease tests have been found to be nonreactive or negative:

Human Immunodeficiency Virus (HIV)

HIV-1/2 Plus O Antibodies (HIV-1/2 Plus O Ab)
Nucleic Acid Test for HIV-1 RNA (HIV-1 NAT)

Hepatitis B Virus (HBV)

HBV Surface Antigen (HBsAg)
HBV Core Antibody (IgG & IgM) (HBcAb)
Nucleic Acid Test for HBV DNA (HBV NAT)

Hepatitis C Virus (HCV)

HCV Antibody (HCVAb)
Nucleic Acid Test for HCV RNA (HCV NAT)

Syphilis*

Rapid Plasma Reagin (RPR) Screen
T. Pallidum IgG

West Nile Virus (WNV)

Nucleic Acid Test for WNV RNA (WNV NAT)

*A donor whose blood specimen is unsuitable for the non-treponemal screening assay, such as the RPR test, or with a reactive result from the non-treponemal screening assay, is cleared for transplantation use only when the result from the treponemal-specific (confirmatory) assay is nonreactive.

Additional tests for other communicable diseases, such as: Cytomegalovirus, Epstein Barr, HTLV I/II, Toxoplasma gondii, T-Cruzi, and Zika Virus, may have been performed. The results of all additional communicable disease tests have been evaluated by the Medical Director and have been found acceptable according to regulations, standards, and VIVEX Biologics policies and procedures.

WARNINGS

The donor of the allograft is screened and tested for relevant communicable diseases and disease agents, and the tissue is microbiologically tested. PRO3 amniotic membrane is processed using aseptic techniques. The allograft is terminally sterilized by a validated electron beam irradiation process in accordance with ANSI/AAMI/ISO 11137. Although efforts have been made to ensure the safety of the allograft, there is no assurance that this allograft is free from all infectious diseases or microbial contamination. As with all allogeneic materials, PRO3 amniotic membrane may transmit infectious agents; however, this risk is greatly reduced by the use of strict donor screening criteria, laboratory testing, aseptic processing, and terminal electron beam irradiation of final product.

DO NOT FREEZE the allograft by any method.

FOR USE IN ONE PATIENT, ON A SINGLE OCCASION ONLY.

DO NOT RE-STERILIZE the allograft by any method. Exposure of the allograft and packaging to irradiation, steam, ethylene oxide, or other chemical sterilant may render the allograft unfit for use.

DO NOT USE ALLOGRAFT IF EXPIRED.

ADVERSE EVENTS AND REACTIONS

Possible adverse events may include:

- Immunologic response (the possibility that a patient may develop alloantibodies should be considered for any patient who might be a future recipient of allograft tissue or cells)
- Transmission of disease of unknown etiology and transmission of infectious agents including but not limited to: HIV, syphilis, or microbial contaminants
- Infection of soft tissue and/or bone (osteomyelitis)
- Fever

STORAGE

PRO3 amniotic membrane must be stored at ambient temperature (2°C to 30°C). It is the responsibility of the tissue dispensing service, tissue distribution intermediary, and/or end-user clinician to maintain the allograft in appropriate storage conditions prior to further distribution or use and to track expiration dates accordingly. Appropriate inventory control should be maintained so that the allograft with the earlier expiration date is preferentially used and expiration is avoided.

80-129 Rev.07
P01-IFU-0004 Rev H

PRECAUTIONS

PRO3 amniotic membrane is processed and packaged using aseptic techniques and terminally sterilized. The allograft must be handled in an aseptic manner to prevent contamination.

Do not use the allograft if either pouch integrity has been compromised.

Use caution when opening, as the allograft is a semi-transparent membrane.

Once the allograft container seal has been compromised, the allograft shall be transplanted promptly, if appropriate, or otherwise appropriately discarded.

The outermost pouch is not sterile and should not be placed on an operative field.

Do not rehydrate PRO3 amniotic membrane prior to use.

ALLOGRAFT PREPARATION

Step 1: Remove the pouch containing the allograft from the box.

Step 2: Inspect the pouches for any holes, tears, or incomplete seals.

Step 3: Using aseptic technique, peel open the outer peel pouch from the chevron end and present the inner pouch to the operative field, when required.

Step 4: Wait to open the inner pouch until ready to place the allograft. Locate the tear notch on the pouch and tear open using caution, as PRO3 amniotic membrane is a semi-transparent membrane.

Step 5: Grasp the allograft and place it directly on the surgical or wound site.

ALLOGRAFT ORIENTATION

The epithelial side of the allograft is facing upwards when a triangle notch is located on the upper left-hand corner of the graft as shown in Figure 1.

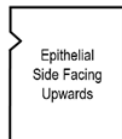


Figure 1.

RECIPIENT INFORMATION

Recipient records must be maintained for the purpose of traceability. It is the responsibility of the tissue dispensing service, tissue distribution intermediary, and/or end-user clinician to provide VIVEX Biologics with information pertaining to the traceability of the allograft used. For this purpose, the postage paid Tissue Utilization Report (TUR) card is provided with the allograft. Once the allograft is used, peel off the small product labels provided on the product packaging and affix on the TUR card and applicable recipient records. Complete the TUR card and mail to VIVEX Biologics, scan, and e-mail to turs@VIVEX.com, or fax to (888) 630-4321.

ADVERSE REACTION OUTCOME AND COMPLAINT REPORTING

Adverse reaction outcomes potentially attributable to the allograft must be promptly reported to VIVEX Biologics at (888) 684-7783. Any other complaints must be promptly reported to Paragon 28, Inc. at (855) 786-2828.

RETURNED GOODS POLICY

Due to the delicate biological nature of a processed allograft, it cannot be returned for credit. If for any reason the allograft must be returned, a return authorization is required from Paragon28, Inc. prior to shipping. It is the responsibility of the healthcare institution returning the allograft to adequately package and label it for return shipment.

The product manufacturer warrants that the allograft will conform to the specifications set forth herein provided that the allograft is handled, stored, and implanted by healthcare providers according to the requirements set forth herein or as provided by it in writing. The product manufacturer makes no other warranties regarding the allograft; specifically disclaims any implied or statutory warranties, including any warranty against disease transmission and infection; and makes no representations or warranties concerning the biological properties or biomechanical properties of the allograft.

Distributed by:



Paragon 28, Inc.
4B Inverness Ct. E, Suite 280
Englewood, Colorado USA
TEL: (855) 786-2828
FAX: (888) 728-1220
www.paragon28.com

Manufactured by:

Vivex Biologics, Inc.
2430 NW 116th Street
Miami, Florida 33167 USA

PRO3-P™ and PRO3-C™ are Trademarks™ of Paragon 28, Inc.



Sterilized by Irradiation